

PART OF THE DURAN® BOTTLE SYSTEM

SWIRLED

DURAN® BAFFLED BOTTLE




**DURAN
WHEATON
KIMBLE**

Excellence in your hands



DURAN® GLS 80® BAFFLED BOTTLE WIDE MOUTH

with GLS 80® thread

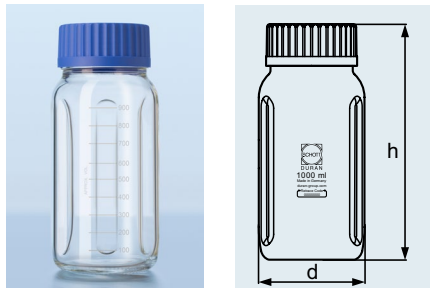

During mixing in standard DURAN® GLS 80® laboratory bottles, the liquid moves in a swirling motion that approximates a solid-body rotation. This is a very inefficient flow pattern, and very little mixing actually occurs. By adding three vertical baffles into the inner surface of the DURAN® GLS 80® bottles, the swirling motion is disrupted and an improvement of the top-to-bottom circulation occurs which produces a greater radial and more effective mixing.

Cat. No.	Capacity (mL)	Thread	d (OD) (mm)	h (mm)	Pack Unit
with screw cap and pouring ring from PP (blue)					
21 286 36 58	250	80	95	110	1
21 286 44 57	500	80	101	153	1
21 286 54 53	1 000	80	101	223	1



> Find your nearest **distributor** on our global network:
www.DWK-LifeSciences.com/DURAN/distributors

