


DURAN® GLS 80® Connection Cap System

screw cap GLS 80®, with four ports,  
GL 18 thread, EPDM seal

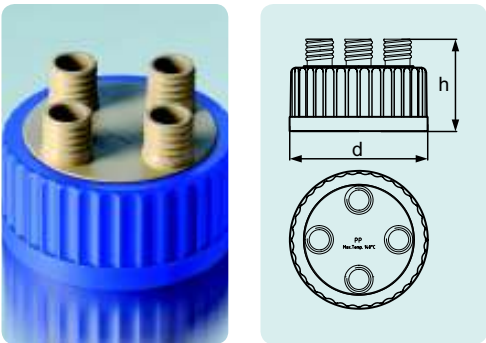


A  
121 °C

Tmax.  
140 °C

Materials used: PP and PTFE. Flexible modular system. Five different tubing diameters (3.2mm; 6.0mm; 8.0mm; 10.0mm and 12.0mm) can be connected. Sterile pressure equalisation is possible by using the syringe membrane filter. Unused ports can be closed with a closed top GL 18 PBT cap.

Typical applications: safe transfer of liquid media within a closed and sterile system (evaporation is reduced).

Cat. No.	Description	Thread	d (OD) (mm)	h (mm)	Pack Unit
293109103	GLS 80® 4-port EPDM seal	80	87	59	2
Accessories					
1160169	Insert for screw cap GL 18, ID 3.2 mm				1
1160170	Insert for screw cap GL 18, ID 6.0 mm				1
1160171	Insert for screw cap GL 18, ID 8.0 mm				1
1160172	Insert for screw cap GL 18, ID 10.0 mm				1
1160173	Insert for screw cap GL 18, ID 12.0 mm				1
1160174	Screw cap for tube connection, blue, GL 18				2
1160167	Pressure equalisation cap set for 4-port, GL 18				1
1170682	Blanking screw cap, red, PBT, GL 18				2

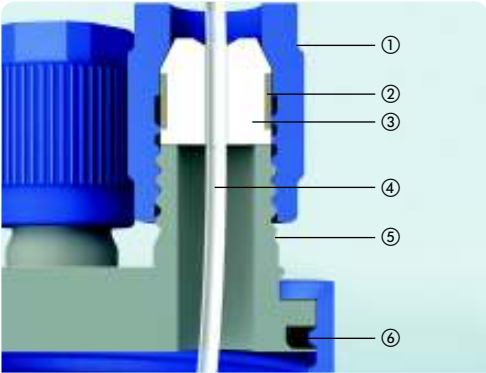

Schematic diagram of GLS 80® connection system

- ① Screw cap GL 18 (PP)
- ② Silicone sealing ring on insert
- ③ PTFE insert / tubing connector
- ④ Tubing (not supplied)
- ⑤ Port (PP)
- ⑥ O-ring seal