


ANALOG MODELS

Stirring Assembly:

Electronic speed control provides constant torque over a speed range from 100 to 1200 RPM, regardless viscosity changes of the medium.

Application:

Stirring viscous gels • Stirring latex mixtures • Mixing polymers, starches & Molasses • Standard reagent mixing • Stir up to 4 Liters of aqueous solution • Aluminium top plate remains cool while stirring • For safe operation, all models include a power indicator light that illuminates when power is applied to the control • Quiet stirring at a gentle 100rpm as well as at a vigorous 1000 rpm operation • Strong magnetic coupling ensures that stir bar remains coupled with drive magnet, even in viscous solutions.

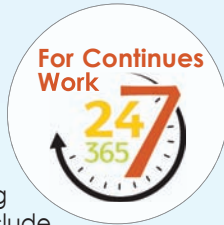

F-23, Magnetic Stirrer, Up To 8 Liter (H₂O)

Powerful drive magnet • Strong motor-stirs up to 8 Liter water • Compact design maximizes use of available bench space • Variable speed control up to 1200rpm • Electronic speed control with "OFF" position • Available as 4-position stirrer, & 6-position stirrer • Powerful motor & strong magnet provide exceptional stirring even with a solution that has a viscosity similar to cooking oil • Durable, easy-to-clean, corrosion-resistant top plate available.



F-23
130x145mm
Stirrer
up to 8Liter

F-33, Magnetic Stirrer, Up To 8 Liter


Similar to F-23 but with bigger plate 18.5 X 17.5 cm.



F-43
305x285mm
Stirrer
up to 30Liter

F-33S, Medium-High Volumes Stirrer, Up To 20Liter

Magnetic stirrer 185 x 175 mm with a heavy duty brushless (non sparking) induction motor with condenser. Exclusive powerful magnet is especially designed for most efficient magnetic coupling. Stir up to 20 Liter water • Gentle or vigorous stirring • Integral ring-stand holder • Include rear bracket for support rod.

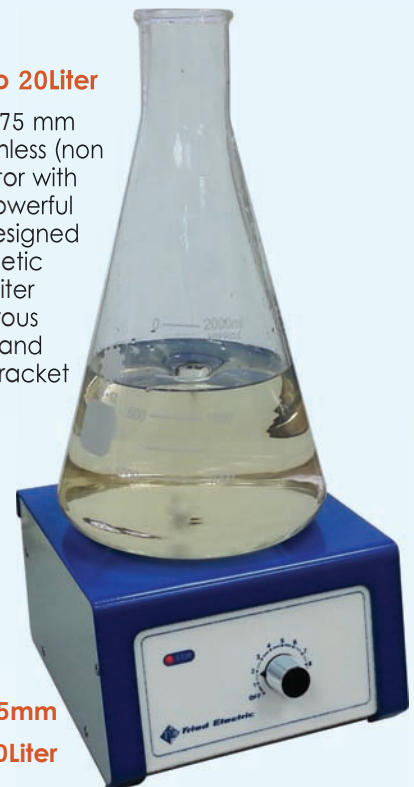

F-33S
185x175mm
Stirrer
up to 20Liter

F-43, Large Volume Stirrer, Up To 30Liter

The giant stirrer is similar to the above described Model F-33 but with a 30.5 x 28.5 cm. top plate. The giant size top plate easily accommodates stirring loads as large as 30 Liter water.

Optional Rod.

Model	Top plate	Capacity	Speed range	Dimensions (WxHxD)	Weight
F-23	W130xD145 mm	8 Liter	100 to 1200 rpm	130x125x186 mm	1.2 kg
F-33	W185xD175 mm	8 Liter	100 to 1200 rpm	130x125x186 mm	1.2 kg
F-33S	W185xD175 mm	20 Liter	100 to 1400 rpm	186x125x197 mm	3.2 kg
F-43	W305xD285 mm	30 Liter	100 to 1400 rpm	305x125x305 mm	4.2 kg